LIFE SOURCE[®]

an **AND** brand
A&D Medical

Wrist Auto-Inflation Blood Pressure Monitor

Trilingual Instruction Manual - Model UB-351
English • Français • Español

AND
A&D Medical

IMPORTANT INFORMATION

Please read this important information before using your UB-351 monitor.

-  Please remember that only a medical practitioner is qualified to interpret your blood pressure measurements. Use of this device should not replace regular medical examinations.
-  Have your physician review your procedure for using this device. He or she will want to verify blood pressure measurements before prescribing or adjusting medications.
-  Consult your physician if you have any doubt about your readings. Should a mechanical problem occur, contact A&D Medical.
-  Do not attempt to service, calibrate, or repair this device.
-  Because your UB-351 monitor contains delicate, high-precision parts, avoid exposing it to extremes in temperature or humidity or to direct sunlight, shock and dust. A&D Medical guarantees the accuracy of this device only when it is stored and used within the temperature and humidity ranges noted on page E-14.
-  Clean the monitor and cuff with a dry, soft cloth or a cloth dampened with water and a mild detergent. Never use alcohol, benzene, thinner or other harsh chemicals to clean monitor or cuff.
-  Remove and replace batteries if monitor is not used for more than six months.

PRECAUTIONS FOR USE

The UB-351 is designed to be used at home, by those who are 18 years and older to monitor blood pressure: systolic and diastolic, and pulse rate.

LifeSource Health Line (Toll-Free): 1-888-726-9966

TABLE OF CONTENTS

WHAT DISPLAY SYMBOLS MEAN....	E-2 - E-3
MONITOR COMPONENTS	E-3
HOW UB-351 WORKS.....	E-4

BEFORE YOU START	E-4
------------------------	-----

TAKING YOUR BLOOD PRESSURE ..	E-5 - E-6
WHAT IS AN IRREGULAR HEARTBEAT	E-7
ABOUT AVERAGE READING.....	E-7
ABOUT PRESSURE RATING INDICATOR™	E-8
ABOUT MEMORY.....	E-9

ABOUT BLOOD PRESSURE	E-10 - E-12
----------------------------	-------------

TROUBLESHOOTING	E-13
SPECIFICATIONS	E-14- E-15

CONTACT INFORMATION	E-15
FRENCH INSTRUCTIONS.....	F-17
SPANISH INSTRUCTIONS	S-35

WHAT DISPLAY SYMBOLS MEAN


Display Symbol	Condition/Error	Recommended Action
	Appears when measuring is in progress, and flashes when pulse is detected.	Remain still and refrain from talking during monitoring process.
	An irregular heartbeat or body movement may have occurred.	Take measurement again and consult with your physician.
 BATTERY FULL	The battery power indicator during measurement.	_____
 BATTERY LOW	The batteries are low when the indicator blinks.	Replace batteries with two new batteries.
	Appears if the pressure value is unstable due to movement during the measurement.	Take another reading, making sure you remain quiet and perfectly still.
	Appears when the pressure reading is unstable due to movement during the measurement.	Take another reading, making sure you remain quiet and perfectly still.
	Appears when inflation is not sufficient.	If your systolic reading exceeds 150 mmHg, then refer to page E-10.
	Appears when the cuff is not fastened correctly.	Check cuff to ensure fit and fastening is correct, and measure again.
 PULSE DISPLAY ERROR	Appears when pulse is not sufficiently detected.	Check for correct fit and position of cuff and take another reading.

LifeSource Health Line (Toll-Free): 1-888-726-9966

WHAT DISPLAY SYMBOLS MEAN

Display Symbol	Condition/Error	Recommended Action
M	Indicates previous measurement stored in memory.	
	Cuff Inflation Meter	Measurement is in progress, remain as still as possible
	Pressure Rating Indicator™	Refer to the section "About Pressure Rating Indicator" for further explanation.

MONITOR COMPONENTS

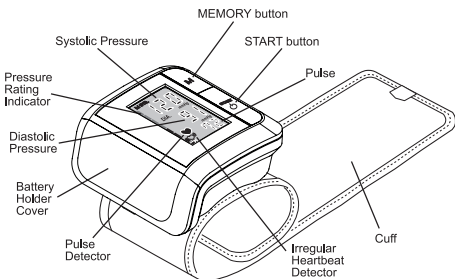In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

HOW UB-351 WORKS

LifeSource® Blood Pressure Monitors are easy to use, accurate and display clear digital measurements. Our technology is based on the “oscillometric method”. The term “oscillation” refers to the measure of vibrations caused by the arterial pulse and is a non-invasive method to determine blood pressure.


Our monitors examine the pulsatile pressure generated by the arterial wall as it expands and contracts with each heartbeat against the monitor’s cuff. This monitor measures your blood pressure while the cuff is inflating.

BEFORE YOU START

The UB-351 provides accurate readings on wrists with a circumference between 5.3" and 8.5" (13.5 cm and 21.5 cm).

The unit requires two "AAA" (1.5 volt) batteries to operate (alkaline batteries recommended). To install batteries (or to replace when "Low Battery" symbol appears), proceed as follows:

1. Release the battery cover as shown to remove from the device.
2. Insert batteries into designated areas ensuring a correct match with positive (+) and negative (-) terminals.
3. Return battery cover by sliding back on to monitoring device.

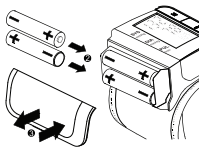NOTE: *Rechargeable batteries will not work with this monitor.*

LifeSource Health Line (Toll-Free): 1-888-726-9966

TAKING YOUR BLOOD PRESSURE

Tips for Blood Pressure Monitoring:

- ☞ Relax for about 5 minutes before measurement.
- ☞ Do not smoke or ingest caffeine at least 30 minutes prior to measurement.
- ☞ Remove constricting clothing and place cuff on bare wrist.
- ☞ Unless your physician recommends otherwise, use left arm to measure pressure.
- ☞ Do not talk during measurement.
- ☞ Do not cross legs. Keep your feet flat on floor during measurement.

Now you are ready. Follow these simple steps:

1. Sit comfortably in a chair near a flat surface.
2. With palm facing upward, attach monitor to the inside of your left wrist and fasten securely. The top edge of the cuff should align with the crease of your wrist.
3. Position your wrist so it is level with your heart, using an arm rest if needed. Positioning of the wrist at heart level is critical for accurate measurement.
4. Press the START button. As the cuff pressurizes, the Cuff Inflation Meter will show on the display screen. It is normal for the cuff to feel very tight.

NOTE: To stop inflation at any time, press START

5. Remain still and refrain from talking or moving during reading.

NOTE: If a reading is not correctly detected, the device will start to automatically re-inflate.

In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

E-5

TAKING YOUR BLOOD PRESSURE

6. After the cuff is fully inflated, it will automatically start to deflate and  symbol will appear to indicate that measurement is in progress. Once the pulse is detected, the symbol will flash with each pulse beat.
7. When measurement is complete, the systolic and diastolic pressure and pulse rate readings are displayed. The cuff will automatically deflate completely, releasing any remaining air.
8. If memory button is pressed when the measurement result is displayed, the reading is not memorized.
9. Press the START button once to turn off power. The device has an automatic power-down feature which turns off power automatically one minute after measurement is obtained.

For the most useful blood pressure history, we recommend that you record the date and time after each measurement. The most accurate history relies not on a single or sporadic readings, but on a consistent monitoring pattern over time.

NOTE: If you have poor blood circulation in winter (especially when your hands and wrists are cold), you may not get a blood pressure measurement.

LifeSource Health Line (Toll-Free): 1-888-726-9966


WHAT IS AN IRREGULAR HEARTBEAT

The UB-351 Blood Pressure Monitor provides a blood pressure and pulse rate measurement even when an irregular heartbeat occurs. The Irregular Heartbeat symbol «♥» will appear in the display window in the event an Irregular Heartbeat has occurred during measurement. An irregular heartbeat is defined as a heartbeat that varies by 25% from the average of all heartbeat intervals during the blood pressure measurement. It is important that you relax, remain still and refrain from talking during measurements.

NOTE: We recommend contacting your physician if you see this «♥» symbol frequently.

ABOUT AVERAGE READING

The UB-351 provides an average blood pressure reading based on the total measurements stored in memory. When taking a measurement, the average blood pressure reading will appear and the number of measurements stored in memory will flash several times before the cuff begins to inflate. The average blood pressure reading will also appear briefly when you recall the measurements stored in memory.

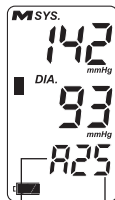
In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

ABOUT BLOOD PRESSURE RATING INDICATOR™

The Pressure Rating Indicator™ is a feature which provides a snapshot of your blood pressure classification based on your measurements. This will let you quickly know what your blood pressure readings mean. Each segment of the bar indicator corresponds to the Seventh Report of the Joint National Committee (JNC7) on Prevention, Evaluation and Treatment of High Blood Pressure from the National Heart, Lung and Blood Institute - May 2003. For a more detailed look at this blood pressure classification, please refer to "Assessing High Blood Pressure" on page E-10.

Pressure Rating Indicator™

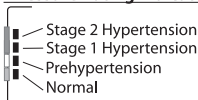Example

Stage 2
Hypertension

The indicator displays a segment, based on the current measurement, corresponding to the JNC7 Classification.

NOTE: Due to other risk factors (e.g. diabetes, obesity, smoking, etc.) in addition to your blood pressure measurement, the Pressure Rating Indicator is approximate. Please consult with your physician for interpretation and diagnosis of your blood pressure measurements.

NOTE: Residents outside of the United States (e.g., Canada and Mexico) should refer to the WHO Classification Table on page E-11 for assessment of their blood pressure measurement.

LifeSource Health Line (Toll-Free): 1-888-726-9966


ABOUT MEMORY

The UB-351 automatically stores up to (90) blood pressure and pulse measurements in its memory. Measurements stored in memory are assigned an index number in the order of the newest to the oldest. The oldest reading displays as "n01".

To retrieve readings, follow these simple steps:

1. When pressing the MEMORY button, the UB-351 will go into Memory display mode.
2. The display will show you the average memory of the readings stored in memory.
3. When you press the MEMORY button again, the total number of memorized data points is shown on the middle line. The UB-351 shows the latest data index number and measurement results sequentially.
4. If the MEMORY button is pressed again, the UB-351 will display the data starting with the chronologically latest data in memory.
5. To retrieve other measurements in memory, press and release the MEMORY button repeatedly until the desired index number is reached (e.g. no.5). The measurement will then be displayed.
6. The display will shut off automatically after two cycles.

If there are no measurements stored in memory, you will see this screen displayed:

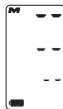To **clear** the measurement history from memory, remove the batteries from the unit and then place the batteries back into the battery holder.

NOTE: If the batteries are removed, the stored data is deleted.

In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

E-9

ABOUT BLOOD PRESSURE

■ What Is Blood Pressure?

Blood pressure is the force exerted by blood against the walls of the arteries. Systolic pressure occurs when the heart contracts; diastolic pressure occurs when the heart expands. Blood pressure is measured in millimeters of mercury (mmHg).

■ Assessing High Blood Pressure for Adults

The following standards for assessing high blood pressure (without regard to age or gender) have been established as a guideline. Please note that other risk factors (e.g. diabetes, obesity, smoking, etc.) need to be taken into consideration and may affect these figures. Consult with your physician for an accurate assessment.

JNC7 Classification Table – for adults within the U.S

BP Classification	Systolic (mmHg)		Diastolic (mmHg)
Normal	less than 120	and	less than 80
Prehypertension	120-139	or	80-89
Stage 1 Hypertension	140-159	or	90-99
Stage 2 Hypertension	greater than or equal to 160	or	greater than or equal to 100

SOURCE: The Seventh Report of the Joint National Committee on Prevention, Evaluation and Treatment of High Blood Pressure for adults. National Heart, Lung and Blood Institute - May 2003.

■ Measurements for those whose systolic pressure is over 150 mmHg

1. Fasten the cuff around your wrist. Raise your arm so that the cuff is at the same level as your heart.
2. Press and hold the START button until a number about 30 to 40 mmHg higher than your expected systolic pressure appears.
3. When the desired number is reached, release the START button. Measurement will begin.

LifeSource Health Line (Toll-Free): 1-888-726-9966

ABOUT BLOOD PRESSURE

WHO Classification Table – for adults outside of the U.S. (e.g. Canada, Mexico)

BP Classification	Systolic (mmHg)		Diastolic (mmHg)
Optimal	less than 120	and	less than 80
Normal	less than 130	or	less than 85
High-Normal	130-139	or	85-89
Stage 1 Hypertension	140-159	or	90-99
Stage 2 Hypertension	160-179	or	100-109
Stage 3 Hypertension	greater than 180 or equal to	or	greater than 110 or equal to

SOURCE: Standards to assess high blood pressure, without regard to age or gender, have been established by the World Health Organization (WHO).

■ Keys to Successful Monitoring:

Blood pressure fluctuates throughout the day. We recommend that you are consistent in your daily measurement routine:

- Measure at the same time every day.
- Sit in the same chair/position.
- Do not cross legs and keep your feet flat on the floor.
- Relax for 5 minutes before measurement.
- Use the correct cuff size to get an accurate reading.
- Sit still during measurement – no talking, eating or sudden movements.
- Record your measurement in a logbook.

In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

E-11

ABOUT BLOOD PRESSURE

■ Establishing Baseline Measurements

The most important method to get an accurate blood pressure measurement is consistency. To get the most benefit out of your monitor, it is important to establish a “baseline measurement.” This helps build a foundation of measurements that you can use to compare against future readings. To build this baseline measurement, devote two weeks for consistent blood pressure monitoring. This involves doing everything the same way when you measure (e.g. measuring during the same time of day, in the same location, sitting in the same chair, etc.). Once you establish your baseline measurement, you can start evaluating if your measurement has been affected based on things like lifestyle changes or medication treatment.

■ How Do I Record My Blood Pressure?

Blood pressure readings are typically recorded with the systolic pressure written first, followed by a slash mark and the diastolic pressure. For example, 120 mmHg systolic and 80 mmHg diastolic measurements are written as 120/80. Pulse is simply written with the letter “P” followed by the pulse rate—P 72, for example. Please see the back of the manual for the blood pressure tracking record.

LifeSource Health Line (Toll-Free): 1-888-726-9966

TROUBLESHOOTING

Problem	Probable Cause	Corrective Action
Nothing appears in the display, even when the power is turned on.	Batteries are drained.	Replace all batteries with new ones. Alkaline batteries recommended.
	Battery terminals are not in the correct position. Battery voltage is too low.	Reinstall the batteries with negative and positive terminals matching those indicated in the battery compartment.
The cuff does not inflate.	Low battery symbol blinks. [If the batteries are drained completely, the mark does not appear.]	Replace all batteries with new ones. Alkaline batteries recommended.
The unit does not measure. Readings are too high or too low.	The cuff is not fastened properly.	Fasten the cuff correctly.
	You moved your arm or body during the measurement.	Make sure you remain very still and quiet during the measurement.
	The cuff position is not correct.	Sit comfortably and still. Make sure the cuff is at the same level as your heart.

NOTE: If the actions described above do not solve the problem, call 1-888-726-9966. Do not attempt to repair the device yourself.

In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

E-13

SPECIFICATIONS

Model	UB-351
Type	Oscillometric
Display	Digital, 10/7-mm character height, Pressure/Pulse displayed simultaneously
Memory	90 readings
Measurement range	Pressure: 20 mmHg to 280 mmHg Pulse: 40 pulses to 180 pulses/minute
Accuracy	Pressure: ± 3 mmHg or $\pm 2\%$, whichever is greater, Pulse: $\pm 5\%$
Pressurization	Automatic, using micropump
Depressurization	Constant air release-valve system
Power source	2 type AAA (1.5 volt) alkaline batteries (included)
Battery life	Approximately 250 measurements
Recommended Operating	
environment	50°F to 104°F (10°C to 40°C) Less than 95% relative humidity
Recommended Storage	
environment	14°F to 140°F (-10°C to 60°C) Less than 95% relative humidity
Wrist Circumference Range	5.3" – 8.5" (13.5 cm – 21.5 cm)
Weight	3.6 oz. (102g) without batteries
Dimensions	Length 1.4" (35mm) Width 2.5" (64mm) Height 2.7" (68mm)

LifeSource Health Line (Toll-Free): 1-888-726-9966

SPECIFICATIONS

Blood pressure measurements determined by the UB-351 are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute for electronic or automated sphygmomanometers.

This LifeSource® blood pressure product is covered by a 5 year Warranty.

CONTACT INFORMATION

**Visit our web site at www.LifeSourceOnline.com
for warranty registration.**

For more information regarding use, care or servicing of your blood pressure monitor, contact:

In USA:

A&D Medical
A division of A&D Engineering, Inc.
1756 Automation Parkway
San Jose, CA 95131
LifeSource Health Line (Toll-Free):
1-888-726-9966
www.LifeSourceOnline.com

In Canada:

Auto Control Medical
6695 Millcreek Drive, Unit 5
Mississauga, Ontario
L5N 5R8 Canada
(Toll-Free): 1-800-461-0991
www.autocontrol.com

In Canada - Auto Control Medical (Toll-Free): 1-800-461-0991

IMPORTANT!

If You Need Assistance with the Set-Up or Operation

We Can Help!

Please call us FIRST before contacting your retailer at

LifeSource Health Line

1-888-726-9966 - Toll Free

Auto Control Medical

1-800-461-0991 - In Canada

A specially trained representative will assist you.

1756 Automation Parkway

San Jose, CA 95131

1-888-726-9966

www.LifeSourceOnline.com

Pressure Rating Indicator™ is a registered trademark of A&D Medical and is patent pending.

© 2011 A&D Medical. All rights reserved.

Specifications subject to change without notice.

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